




Complete Vacuum Delivery System

All-In-One, One For All



Kiwi's flexible stem and low profile cup facilitates proper placement & deliveries with **ALL** fetal head positions.

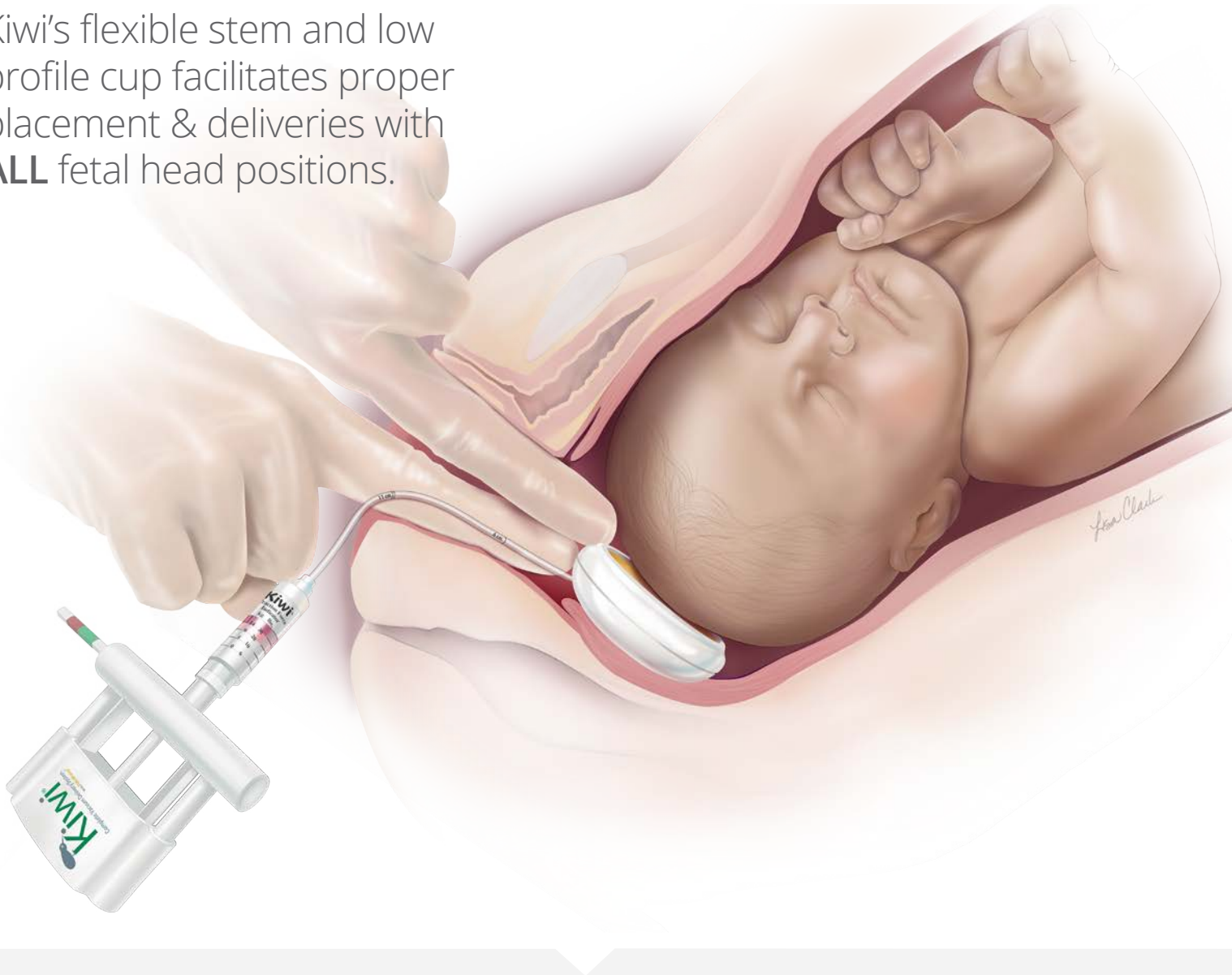

All Positions, All Encompassing, All-In-One

omnicup: | amni-'kəp |
1. a vacuum delivery device that can be used in **all** fetal head positions.

All Fetal Head Positions

The OmniCup's maneuverable design allows it to be used regardless of fetal head position.

All Encompassing


The OmniCup handles any vacuum delivery situation.

All Components in One

The OmniCup is self-contained, eliminating any assembly & other components.



Align Smallest Diameter of Head

Unless the fetal head is presenting directly OA, a flexible stem is needed to position the cup over the optimal location (flexion point). This is necessary in order to flex the baby's head, which presents the smallest diameter of the head in line with the birth canal.

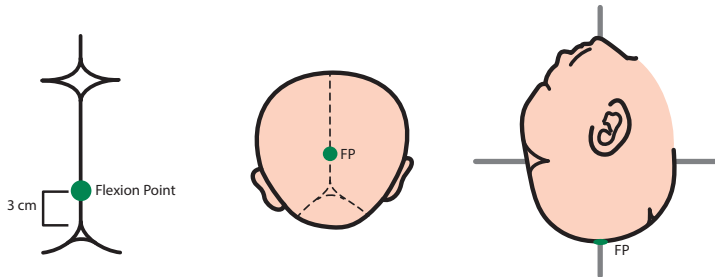

Proper Placement With Flexible Stem

The key to any vacuum delivery is proper cup placement. Kiwi's low profile cup and flexible stem help ensure the cup is placed properly: allowing placement over the flexion point in **ALL** fetal head positions.

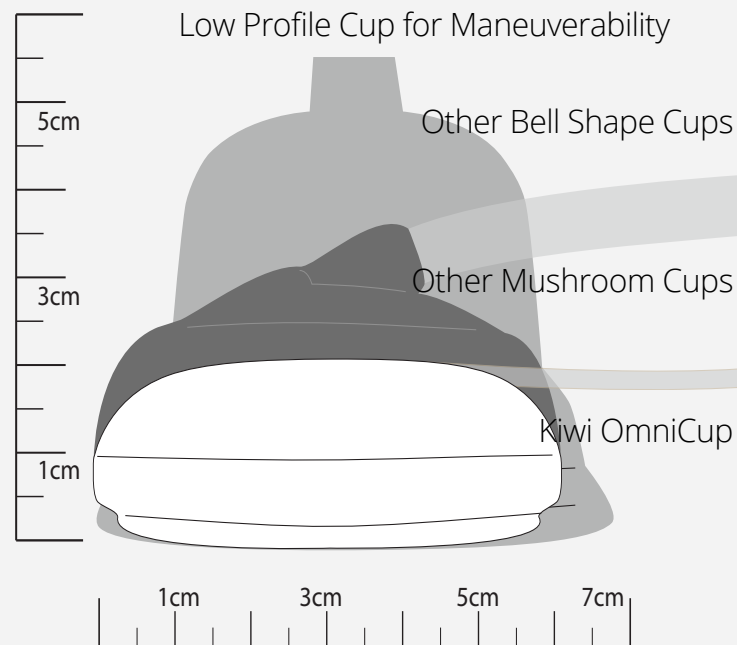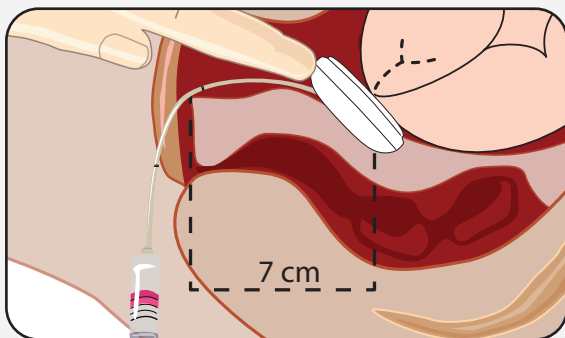

Product Support

With Kiwi you will get more than a cup! Clinical Innovations provides a **5-Step Vacuum Technique** training to improve vacuum assisted delivery. Ask about a hands on simulation workshop.

Ordering Information

Product Description	Product Code	Box Quantity	Positions
Kiwi OmniCup with Traction Force Indicator	CI VAC-6000MT	5/ea	OA, OP, OT, C/S
Kiwi OmniCup	CI VAC-6000M	5/ea	OA, OP, OT, C/S
Kiwi Omni-C Cup	CI VAC-6000C	5/ea	C/S
Kiwi ProCup	CI VAC-6000S	5/ea	OA

OA= Occiput Anterior OP= Occiput Posterior OT=Occiput Transverse C/S=Cesarean Section

Local Distributor:

clinicalinnovations.com.au

36 Matheson Street, Virginia
Queensland 4014 | Australia
Phone +07 3862 2000 | Fax +07 3862 2077
enquiries@clinicalinnovations.com.au

